

MEDIKA[®]

SILKROLL

SPECYFICATION

SILK FEEL



The device with an internal roller is a non-invasive mechanical treatment that utilizes compression and microvibrations. Its operating principle is based on silicone balls rotating around the roller in a 360° range, generating compression and microvibrations. During the rotation, the balls exert pressure on the skin, creating a „pulsating compression“



Maximum
slimming

effect. As a result, there is a continuous kneading and pushing movement in both directions, and the tissue experiences pressure and a lifting effect without causing compression or damage to the skin.

The device features two rollers:

- ✔ The larger roller, with 50 silicone balls, is designed for body care.
- ✔ The smaller roller, with 60 silicone balls, is intended for facial and neck care.

The silicone/metal balls are arranged in a honeycomb pattern with specific density and diameter. The rotation direction and pressure ensure effective transmission of micro-compression to the tissues. The unique handle generates mechanical compression and microvibrations on the skin, which cannot be replicated manually or by other techniques.


Compression and microvibration therapy replaces the traditional „suction-stretch massage” method, using micro-compression instead of suction and microvibrations instead of stretching.



Mode of Operation


Pressure on the tissues stretches the cells, naturally and deeply stimulating their activity, blood flow, and oxygenation. Fat deposits subjected to pressure are loosened, leading to their breakdown and a reduction in cellulite. Additionally, pressure on deeper muscle groups fully softens and stretches them, reducing muscle stiffness, relieving pain, and accelerating metabolism.

The device eliminates fluid retention and accumulations in the tissues, improving their condition and re-firming the skin, effectively shaping the body. Additionally, it stimulates fibroblasts, increasing collagen and elastin production, while also enhancing blood flow and oxygen delivery. As a result, wrinkles are smoothed, puffiness and under-eye bags are reduced, and the skin appears rejuvenated and firmer. This technology has been clinically proven effective in toning muscles and sculpting the body and face, promoting rejuvenation, as well as enhancing shape and firmness, including the chest area.

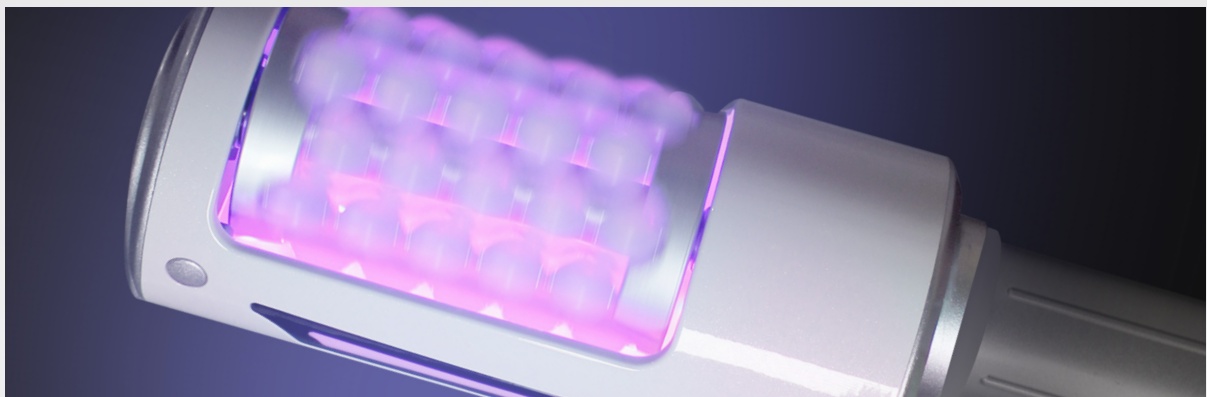


Sculpting effect

It shapes the body and smooths wrinkles.

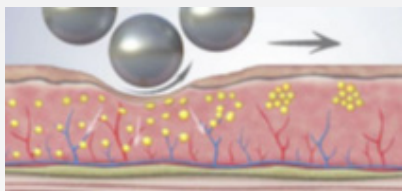


Effect on Muscles

It tones the muscles and firms the tissues.

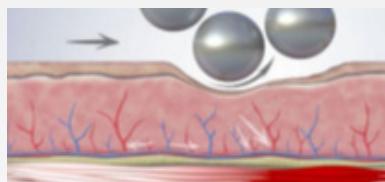


Vascular Effect:

It improves blood circulation and nourishes the tissues.

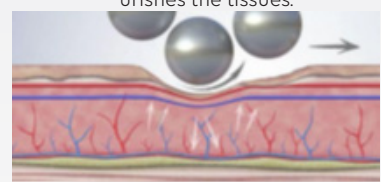


Soothing Effect:

It enhances oxygenation and alleviates

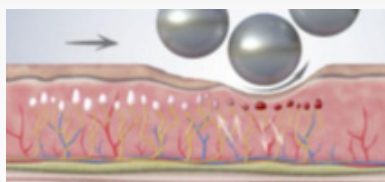

Draining Effect:

It eliminates toxins and supports lymphatic drainage.

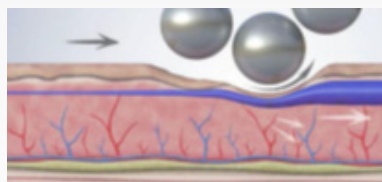



Application

- Improves lymphatic drainage and circulation
 - Supports deep detoxification
 - Relieves muscle stiffness and pain
- Breaks down subcutaneous fat deposits
Eliminates cellulite and shapes the body

Although the effects are visible after the first treatment, a series of sessions (usually 6-12) is typically recommended to achieve optimal results. The duration and frequency of the sessions will depend on the individual needs and goals of the patient. After several treatments, the body can reduce its volume by 1-2 sizes!



Face:

- Improves the formation of blood vessels
 - Supports the production of collagen and elastin
- Breaks down subcutaneous fat deposits
— Eliminates cellulite and shapes the body

— Relaxes muscles, reduces expression wrinkles, and regulates muscle tension.

As a result, we achieve:

- Reduction of puffiness
- Combatting tissue sagging
- Visible rejuvenation and firming of the facial skin



Device Properties

The unique 360° rotating head with an intelligent drum allows for continuous, long-lasting operation in a safe and stable manner.

The head features an LED display showing the time and speed, along with a backlit LED column, which facilitates the control and adjustment of the rotation direction and speed on the body handle.

One-button switching between forward and reverse rotations.

The silicone balls are **flexible and smooth**, providing a gentle rolling process that does not cause discomfort. The movement is soft

evenly distributed, allowing for an effective massage and lifting to achieve the best results. There is no need for a time-consuming massage by a cosmetologist – **the operation is simple and safe.**

THREE COLORS TO CHOOSE FROM.

The heads are equipped with LED lights in red, blue, and purple.

RED LIGHT

It has anti-aging effects, stimulates collagen and elastin production, reduces wrinkles, increases firmness, and improves skin elasticity. It enhances circulation, revitalizes metabolic processes, and stimulates the regeneration of both deep tissues and the skin.

BLUE LIGHT

It helps to calm the sebaceous glands, and the reduction in their activity leads to fewer breakouts. It soothes sensitive skin with enlarged capillaries. It also stimulates the body to increase its activity.

PURPLE LIGHT

In other words, the combination of red and blue light – thanks to the combination, the light becomes dual-spectrum, with antibacterial properties and stimulates melanin production.





Technical Specifications

WEIGHT

20 kg

DIMENSIONS

97x40x40 cm

SCREEN

touchscreen 8"

NUMBER OF HEADS [PCS]

2

POWER CONSUMPTION

50 VA

ROTATION SPEED (SMALL)

Level I — 200 Level IX - 790

(LARGE) ROLLER

Level I — 145 Level IX — 580

SYSTEM

Light Therapy

POWER SUPPLY

230 V AC, 50Hz

DEGREE OF PROTECTION

20 IP

WARRANTY

30 months

MEDIKA[®]

Potokowa 47 Street
54-105 Wrocław

kontakt@medikapoland.pl
71 350-47-65 | 794-151-415